



LAB LINE

Stereo zoom microscope with or without halogen illumination, for the laboratory, training centres, quality control or agriculture

Features

- The KERN OZL-45 stereo zoom microscope series will impress you with its excellent optical characteristics, easy operation and high level of ergonomic working comfort
- The Halogen incident and transmitted illumination included as standard guarantees the very best illumination of your sample
- The high-quality optics, together with a large working surface offers the highest level of comfort for your applications
- The zoom objective offers you continuous magnification from 7,5×–50×
- The KERN OZL-45 series is available as a binocular version
- The pillar stand offers you the highest level of flexibility and the freedom to remove the microscope head and to integrate it into other modular systems, for example into a universal stand

- A large selection of eyepieces, external illumination units as well as auxiliary objectives are available as accessories
- A protective dust cover, eye cups, as well as multi-lingual user instructions are included in the scope of delivery
- Please find detailed information in the following model outfit list

Scope of application

- In vitro fertilisation, detection of parasites, zoology and botany, tissue preparation, section, quality control

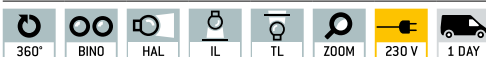
Applications/Samples

- Samples with focus on three-dimensional impression, zoom with variable magnification (depth, thickness), e.g. insects, seeds, circuit boards, components

Technical data

- Optical system: Greenough optics
- Brightness adjustable
- Tube 45° inclined
- Magnification ratio: 6,7:1
- Interpupillary distance 55 – 75 mm
- Diopter adjustment: Both-sided
- Overall dimensions W×D×H 330×270×460 mm
- Net weight approx. 5 kg

STANDARD

























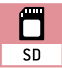
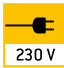



Model	Standard configuration					
	Tube	Eyepiece	Field of view mm	Objective Zoom	Stand	Illumination
KERN						
OZL 45 1	Binocular	HSWF 10×/ø 23 mm	ø 33 – 5	0,75× – 5,0×	Pillar style	12 V/10 W Halogen (incident) 12 V/10 W Halogen (transmitted)

Eyepiece	Specifications - Objectives				
	Magnification	Standard 1,0×	Auxiliary objectives		
			0,5×	0,75×	2,0×
HWF 5×	Total magnification	3,75× - 25×	1,875× - 12,5×	2,813× - 18,75×	7,5× - 50×
	Field of view mm	∅ 31 - 4,6	∅ 61,3 - 9,2	∅ 41,3 - 6,1	∅ 16 - 2,5
HSWF 10×	Total magnification	7,5× - 50×	3,75× - 25×	5,625× - 37,5×	15× - 100×
	Field of view mm	∅ 33 - 5	∅ 65 - 10	∅ 44 - 6,7	∅ 16 - 2,5
HWF 15×	Total magnification	11,25× - 75×	5,625× - 37,5×	8,438× - 56,25×	22,5× - 150×
	Field of view mm	∅ 24 - 4,2	∅ 48 - 8,5	∅ 32 - 5,6	∅ 12 - 2
HSWF 20×	Total magnification	15× - 100×	7,5× - 50×	11,25× - 75×	30× - 200×
	Field of view mm	∅ 20 - 3,5	∅ 40 - 7	∅ 26,7 - 4,7	∅ 10 - 1,8
HWF 25×	Total magnification	18,75× - 125×	9,375× - 62,5×	14,063× - 93,75×	37,5× - 255×
	Field of view mm	∅ 15,8 - 2,4	∅ 31,5 - 4,8	∅ 24,1 - 3,2	∅ 7,9 - 1,2
Working distance		113 mm	177 mm	117 mm	35 mm
Maximum sample height		120 mm	60 mm	90 mm	165 mm

Model outfit		Model KERN	Order number	
		OZL 451		
Eyepieces (30,0 mm)	HWF 5×/∅ 23,2 mm	○ ○	OZB-A4112	
	HSWF 10×/∅ 23 mm	✓ ✓	OZB-A4118	
	HWF 15×/∅ 15 mm	○ ○	OZB-A4119	
	HSWF 20×/∅ 14,5 mm	○ ○	OZB-A4120	
	HWF 25×/∅ 11,7 mm	○ ○	OZB-A4121	
Auxiliary objectives	0,5×	○	OZB-A4509	
	0,75×	○	OZB-A4510	
	2,0×	○	OZB-A4206	
Stand	Pillar style, with 12 V/ 10 W Halogen Illumination (transmitted + incident)	✓		
Stage plate	Frosted glass/∅ 95 mm	✓	OZB-A4805	
	Black-white/∅ 95 mm	✓	OZB-A4806	
Illumination	12 V/ 10 W spare bulb (transmitted + incident)	✓	OZB-A4804	
Mechanical stage (Pre-assembling on request)	Stage size W×D 180×155 mm, Travel 75×55 mm, for transmitted and incident illumination	○	OZB-A4605	
External illumination	Please find the information about external illumination units in the catalogue on page 83 and on our website www.kern-sohn.com			

✓ = Included with delivery

○ = Option

 360°	360° rotatable microscope head	 FL-LED	Fluorescence illumination for compound microscopes With 3 W LED illumination and filter	 WLAN	WLAN data interface: For transmitting of the picture to a mobile display device
 MONO	Monocular Microscope For the inspection with one eye	 PH	Phase contrast unit For a higher contrast	 HDMI	HDMI digital camera For direct transmitting of the picture to a display device
 BINO	Binocular Microscope For the inspection with both eyes	 DF	Darkfield condenser/unit For a higher contrast due to indirect illumination	 SOFTWARE	PC software To transfer the measurements from the device to a PC.
 TRINO	Trinocular Microscope For the inspection with both eyes and the additional option for the connection of a camera	 POLAR	Polarising unit To polarise the light	 AUTO ATC	Automatic temperature compensation For measurements between 10 °C and 30 °C
 ABBE	Abbe Condenser With high numerical aperture for the concentration and the focusing of light	 INFINITY	Infinity system Infinity corrected optical system	 IP	Protection against dust and water splashes IPxx The type of protection is shown by the pictogram.
 HAL	Halogen illumination For pictures bright and rich in contrast	 ZOOM	Zoom magnification For stereomicroscopes	 BATT	Battery operation Ready for battery operation. The battery type is specified for each device.
 LED	LED illumination Cold, energy saving and especially long-life illumination	 PARALLEL	Parallel optical system For stereomicroscopes, enables fatigue-proof working	 RECHARGE	Battery operation rechargeable Prepared for a rechargeable battery operation
 IL	Incident illumination For non-transparent objects	 SCALE	Integrated scale In the eyepiece	 230 V	Mains adapter 230V/50Hz in standard version for EU. On request GB, AUS or USA version.
 TL	Transmitting illumination For transparent objects	 SD	SD card For data storage	 230 V	Power supply Integrated in balance. 230V/50Hz standard EU. More standards e.g. GB, AUS or USA on request.
 FL	Fluorescence illumination For stereomicroscopes	 USB 2.0	USB 2.0 digital camera For direct transmitting of the picture to a PC	 1 DAY	Package shipment The time required to manufacture the product internally is shown in days in the pictogram.
 FL-HBO	Fluorescence illumination for compound microscopes With 100 W mercury lamp and filter	 USB 3.0	USB 3.0 digital camera For direct transmitting of the picture to a PC		

Abbreviations

C-Mount	Adapter for the connection of a camera to a trinocular microscope	LWD	Long Working Distance	SWF	Super Wide Field (Field number at least \varnothing 23 mm for 10 \times eyepiece)
FPS	Frames per second	N.A.	Numerical Aperture	W.D.	Working Distance
H(S)WF	High (Super) Wide Field (Eyepiece with high eye point for wearers of glasses)	SLR Kamera	Single-Lens Reflex camera	WF	Wide Field (Field number up to \varnothing 22 mm for 10 \times eyepiece)

Your KERN specialist dealer: